

# Introduction to Population Receptive Field models

Sam Schwarzkopf

Before the start of the actual cookbook on using SamSrf for population receptive field (pRF) mapping, here is a basic introduction to how our method works. Our main procedure is based on Dumoulin & Wandell, 2008, *NeuroImage*, 39(2):647-60. Other analyses are also possible but are only mentioned in passing in this cookbook. Please refer to the `ModelHelp` tool for more details.

## What is a pRF?

As any textbook on neurophysiology will tell you, a neuron's receptive field is the region of space that can excite the neuron when we put a stimulus in it. "Space" can mean a lot of things in this context. Most frequently, people will mean *visual space*, that is, the retinotopic location. However, there are other stimulus spaces, such as maps of the frequency of a tone or locations on your skin. In fact, a visual receptive field is just a tuning function for position in visual space. So in essence what we want to do is map tuning functions for a particular stimulus space – and this could be anything, be it visual position, sound frequency, selectivity to particular objects, your location in the environment, or even things like language. So, don't let anyone tell you that pRF models are just a fancy form of retinotopic mapping. There is potential for so much more.

Naturally, we are not measuring single neurons but the activity of a large number of neurons combined (and, in fMRI, through an indirect, metabolic proxy measure of neuronal activity). This is why we call it a *population* receptive field. The properties of the pRF are not necessarily an up-scaled version of the receptive fields of single neurons. The pRF will contain those parameters but is also dependent on the variability in the position of the individual receptive fields, how much visual space the receptive fields of all the neurons cover, and it also incorporates the extra-classical or contextual interactions between neurons within (and outside) the population.

Last but not least, the acronym is spelled with a lower-case *p* for some arcane reasons. This is probably in reference to the fact that in the neurophysiology world there are such things as *classical* receptive fields that are abbreviated as cRF. For all I care you can spell it PRF.

## How do we map pRFs?



Simply put, for each population of neurons in the brain we want to estimate the parameters of the pRF from the fMRI (or other?) data we collected. The parameters of the pRF can be many things. In a simple case it would be the position (in Cartesian coordinates,  $x_0$  and  $y_0$ ) and the spatial spread of the pRF (since we use a 2D-Gaussian model, this is called  $\sigma$ , i.e. *sigma*). But more complex receptive field models are also possible, such as antagonistic centre-surround profiles (modelled as a difference-of-Gaussians function), or oriented and/or asymmetric profiles. Putting together all the parameters for neighbouring points on the cortical sheet gives us a map for this parameter.

In order to estimate these parameters, we need to have two sets of information: the observed data (i.e. the fMRI time series) and the time series of the stimulus. In SamSrf the latter is called the **stimulus apertures**, but that's just a name for representation of the stimulus we presented to participants during the scan. In fact, it is usually a simplified version of the stimulus that just encodes whether a stimulus was present at each location in the visual field at a given time. In

earlier SamSrf versions (and to our knowledge in other toolboxes), the aperture movie is used directly as input to pRF models. However, as of SamSrf 9 we need to use **vectorised apertures**. This means that we define a set of stimulus locations (such as the pixels in the aperture movie), each of which is assigned a coordinate in visual space. The aperture matrix contains locations in rows, and time points (fMRI volumes) in columns. The reason for doing it this way is that it makes the design of apertures much more flexible (you can use non-symmetric apertures or even just a line of pixels for one-dimensional tuning models). Coordinates can also be defined directly in stimulus space (e.g. degrees of visual angle) rather than the image/aperture space. Most importantly though, this approach speeds up the analysis considerably – in our hands usually by a factor of 3-5!

Armed with pRF model and the stimulus apertures we can now fit the pRF model. What we need to do is find the pRF parameters (like position and size) that can best explain the observed time series. We can predict the time series a particular receptive field profile would produce if faced with the apertures of our stimulus sequence. At each moment in time (i.e. each TR of our fMRI scan) we calculate the overlap between the stimulus aperture and the receptive field profile. As already described, the apertures are just a representation of the stimulus in space and time. The receptive field profile is simply the predicted response for a given location. We multiply this prediction, pixel-by-pixel, with each time point of the aperture. Subsequently we calculate the mean of this pixel-wise overlap to estimate the neural response. Since SamSrf 7, this takes into account the pRF size but in previous versions it simply calculated the overlap across the whole aperture frame. To our knowledge, other pRF tools do the latter but arguably the former is more biologically correct – practically this isn't a big issue though, because it only affects the beta amplitude but doesn't change your pRF shape. If required, changing this is also very straightforward.

As a toy example, consider two pRFs: You see one is in the upper right visual field (coordinates 0.5, 0.5) and very large, while the other is in the upper left visual field (-0.5, 0.5) and very small. Next you need the stimulus apertures. You will have one if you run a pRF mapping experiment but if you don't, ask one of the people who have done the analysis already. This allows us to predict a neural response which will look something like this:

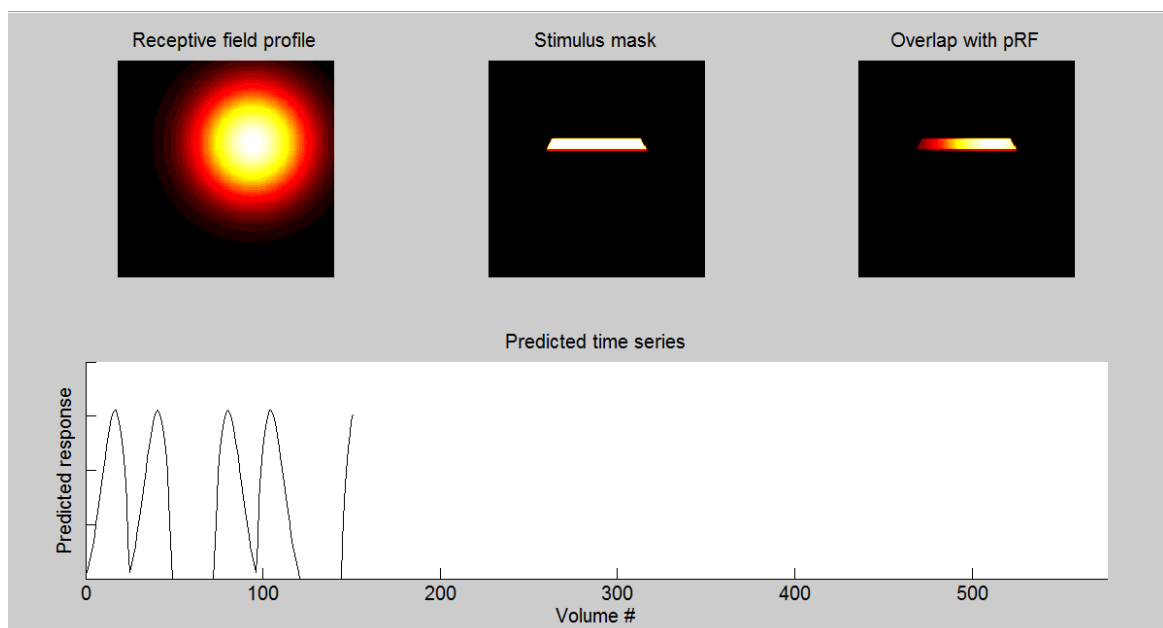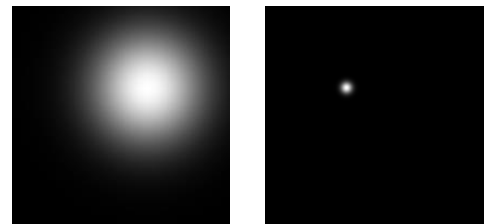

In the top row you see the receptive field profile on the left, the stimulus aperture at the current time point in the middle, and the overlap between the two on the right. In the bottom row we see the predicted response of this hypothetical receptive field. Of course, because we are measuring BOLD signals with fMRI our time series does not reflect the underlying neural activity directly. Therefore, we must convolve the predicted time series with the haemodynamic response function (HRF). For this we can either take a canonical HRF or measure the HRF specifically for each subject.

Theoretically, using a subject-specific HRF is probably better because it takes the individual differences in the HRF into account. In practice, this does not seem to make much of a difference though as Dumoulin & Wandell showed in their seminal 2008 study. Naturally, any individual HRF is only as good as the data that you used to model it. Using a bad HRF can completely screw up your pRF model. Another, more efficient approach is to first fit a pRF based on a canonical HRF, then use the predicted time series as predictor for fitting the individual HRF, and then refitting the pRF model with this individual HRF. An advantage of this approach is that it uses all the pRF mapping data you acquired to fit the HRF rather than spending extra time on measuring the HRF. This approach can also improve the goodness of fit of your pRF models – but this is also somewhat circular. Whether or not this circularity is a problem we don't know, and Sam doesn't have time to find out (if you're interested in doing that study, get in touch and we are happy to consult...). Either way, in many scenarios the canonical HRF is completely acceptable, and most people probably do not need to worry about this.

### Fitting pRF models

The tricky thing we must do now is to find the receptive field profile (or more precisely its associated parameters) that can produce a time series that matches the observed time series. What we therefore need to do is fit a model with several free parameters (in this case at least three) to the observed data. Ideally, we want to do this straight out of the box but this is tedious and difficult. So, what we can do is apply a coarse-to-fine approach in which we first find a way to disambiguate the parameters and then refine them further.

The way this is done for most models in SamSrf is to run an extensive grid search, the *coarse fit*, where we first generate a large search space with plausible combinations of the free parameters. For each point in this multidimensional search space we generate a predicted time series and correlate this with the actually observed data. The search parameters that produced the best correlation are then used in the *fine fit* where we run an optimization procedure to further hone the fitted parameters to achieve the (hopefully) optimal fit. The good thing about the coarse fit is that we can do this simultaneously for many vertices at the same time, which makes the procedure quite fast. (Unfortunately, this appears to only work in MATLAB R2011b and higher. On older versions you must do it at a snails-pace one vertex at a time which takes ages...). As of SamSrf 8, several alternatives for the standard fine fitting are available (see next section for details).

The granularity and dimensions of this search space for the coarse fit you can define when you specify your pRF model. The more predictions this coarse fit contains, the longer it will take to generate predictions and how much disc space and memory the search space requires. There is probably a sweet spot between having a well-sampled search space and unnecessarily wasting too much time and memory on the coarse fit. Obviously, a pRF model with more free parameters also requires more resources and will put limits on how accurately you can sample the search space. **Importantly, you only really need to generate this search space once, provided you collect the**

same amount of data in the same sequence for each subject. If the stimulation paradigm differed between subjects, then unfortunately you also need individual apertures and thus an individual search space.


For both the coarse and fine fit, the best fitting model is determined by the Pearson correlation between the observed and predicted time series. It is quantified by the coefficient of determination (the “goodness of fit”), which is the square of the correlation coefficient. Once the parameter fitting is complete, SamSrf subsequently runs a regression analysis to fit a *beta* parameter to estimate the response amplitude, and a baseline (intercept) parameter to estimate the baseline. Since we typically z-standardize the observed time series, the baseline will typically be close to zero but it can vary. More importantly, the beta can be positive or negative. Negative betas are often a sign of some artifactual data, such as regions outside the mapping stimulus or blood vessel artifacts.

### Algorithms for fine fitting

By default, the fine fit uses an optimisation procedure known as the Nelder-Mead simplex search algorithm, which is implemented in MATLAB’s Optimization Toolbox under the name *fminsearch*. This creates a little polytope around the seed parameters with  $k+1$  vertices, where  $k$  is the number of free parameters (so for a 1D tuning curve with peak and tuning width as parameters, this would be a triangle). It then moves around and shrinks this polytope as required to home in on the local minimum. You can read more about this here:

[https://en.wikipedia.org/wiki/Nelder%E2%80%93Mead\\_method](https://en.wikipedia.org/wiki/Nelder%E2%80%93Mead_method)

This algorithm is very effective, but it can be slow and it also has some problems. It may also be overkill for most pRF applications, at least standard retinotopic mapping purposes. SamSrf includes a function called *samsrf\_gsr2* which allows you to visualise the correlations for each grid position in the search space. Here is an example for a standard 2D Gaussian pRF. Each plot shows the correlations for the  $x_0$  and  $y_0$  positions in the search space for a given pRF size (Sigma):

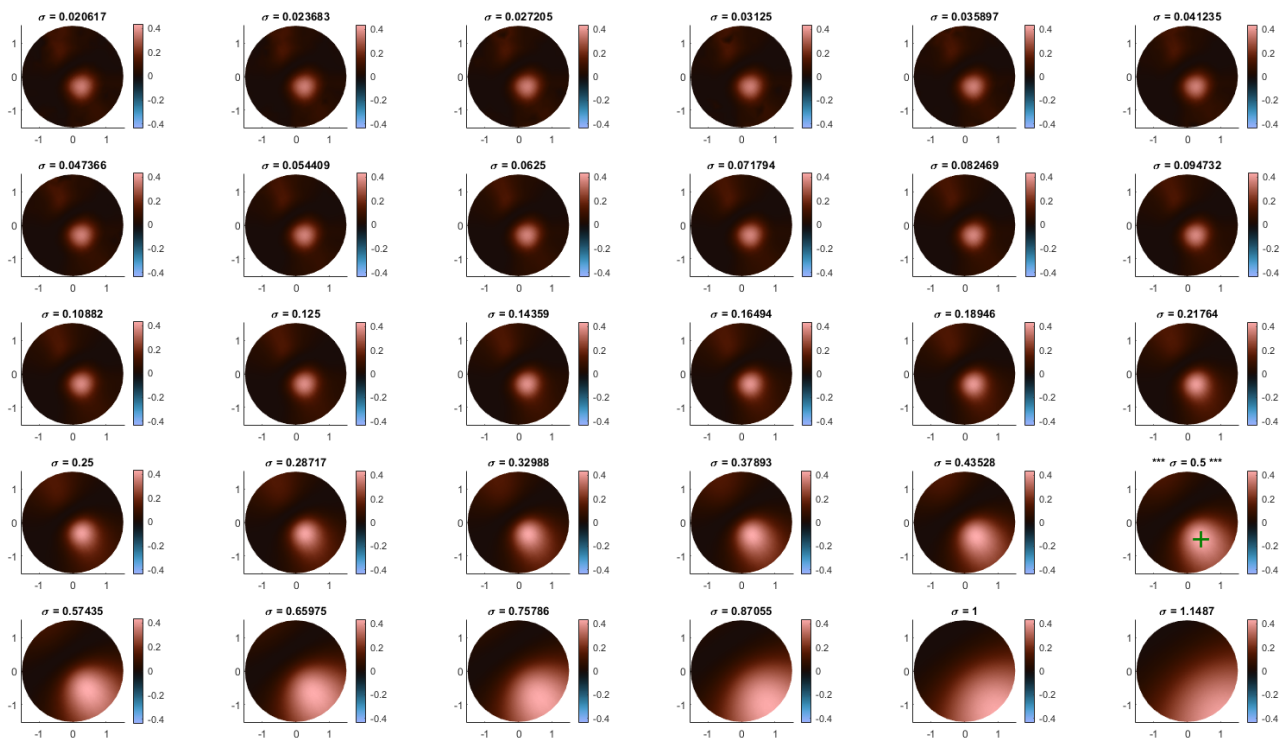

What this shows is that there is a clearly defined cluster of high correlations (pink blob in the bottom right quadrant). The peak correlation is pretty consistent between the various pRF sizes. Only for very large pRFs the cluster becomes blurry and smeared peripherally. The maximum correlation across the whole search space is indicated by the green cross (rightmost column, second to bottom row). This shows that there is actually a very smooth gradient in terms of best fits and very little ambiguity about the parameters. Algorithms like the Nelder-Mead method are really designed to deal with much more complex problems. So for many purposes we can optimise the parameters using much simpler algorithms, which are available in SamSrf:

1. *Slow Coarse Fit*: Instead of only taking the best correlating prediction from the search space it averages the parameters of a top percentile you can define (*Model.Coarse\_Fit\_Percentile*). This can yield pretty good pRF estimates while being considerably faster than the fine fit. But it is slower than the standard coarse fit. There is probably an optimum for what search space granularity is required to obtain sufficiently precise estimates. Simulations can help here. Also keep in mind that this method is necessarily constrained by the limits of your search space. If you expect pRFs to fall outside the stimulus range your search space should extend considerably beyond that and you may need to include  $x_0, y_0 = 0$  (a typical polar search grid would not contain that).

2. *Adjustable parameter tolerance*: You can also adjust the parameter tolerance of the Nelder-Mead algorithm. If the tolerance is higher than default ( $1e-4$ ) this trades accuracy for speed. This was the very problem we ran into with SamSrf 6 because its tolerance was too lenient ( $1e-2$ ) for our purposes. But for many applications this may be completely acceptable & you could also try intermediate tolerance levels. Simulations could help here, too.

3. *Hooke-Jeeves algorithm*: You can also use Hooke-Jeeves pattern-search as an alternative to the Nelder-Mead algorithm. This algorithm works by stepping away from the seed parameters along each dimension. If a better fit is found, the current parameter is moved there. If not, the step size is shrunk by half. See here for more details:

[https://en.wikipedia.org/wiki/Pattern\\_search\\_\(optimization\)](https://en.wikipedia.org/wiki/Pattern_search_(optimization))

This algorithm is comparably fast but not as precise and robust as the Nelder-Mead method. But as we have seen, for a lot of pRF data a simple approach may suffice. For speed, SamSrf only runs a small number of iterations for this fit. (Note: this only applies to pRF forward-modelling; when fitting reverse correlation pRFs/CFs the algorithm uses the default number of iterations). This method seems to perform well for standard 2D pRF models (especially compared to a lenient Nelder-Mead tolerance) but that will probably depend a lot on your situation. Again, this might be worth running some simulations.

## Parallel processing

**SamSrf relies strongly on parallel computing. This means that if you have the Parallel Computing Toolbox in MATLAB installed it will take advantage of multi-core computers and can in theory also be used on high-performance clusters.** This will speed up several computationally expensive procedures, such as smoothing, back-projection, and of course the slow fine-fitting algorithm. SamSrf may be very slow without it. In SamSrf 8, we removed support for non-parallel processing completely, but if you want to use it without it you can in theory edit the functions accordingly.

## Analyses available in SamSrf

Below follows a list of pRF methods currently available in SamSrf. Please also refer to the [ModelHelp](#) tool for more information about specifying models and running these analyses:

1. *Simple Gaussian*: This is the most basic pRF model that is most widely used in the literature. It models the pRF as a two-dimensional Gaussian. It is defined by its position ( $x_0$  and  $y_0$ ) and its size (sigma).
2. *Difference-of-Gaussians*: This procedure models centre-surround antagonistic receptive fields. In addition to position ( $x_0$  and  $y_0$ ), we have two sigmas (one for the centre, one for the surround) and a parameter for the ratio of the amplitudes (delta) of the two components. The second surrounding Gaussian is subtracted from the first, the centre Gaussian. This function produces a Mexican hat shape. **Importantly, you cannot directly translate the centre sigma into a classical sigma from the standard simple Gaussian pRF. Rather, you need to quantify the full-width at half-maximum of this function to determine the pRF size.**
3. *Multivariate Gaussian*: This is an oriented elliptical Gaussian so this has two pRF size parameters (sigma1 and sigma2, for the horizontal and vertical axis, respectively) and an orientation parameter (phi). **There is some ambiguity in this model because the model is rotationally symmetric. You need to account for this in your post-modelling analysis.**
4. *Gaussian Curve*: This procedure is very similar to the standard pRF approach except that it is fitting a one-dimensional tuning curve rather than a two-dimensional receptive field. The apertures in this case should contain bars traversing the stimulus space in only one direction. This direction can be either horizontal or vertical. The tuning curve is modelled as a Gaussian pRF that is fixed in the orthogonal dimension (so if the stimulus moves horizontally, the pRF's vertical position is fixed at 0. **The assumption here is that the stimulus space is linear. So if your stimulus dimension is non-linear you must log-transform it first** (e.g. in tonotopy or spatial frequency).
5. *Polar Angle Curve*: This is the same as the Gaussian Curve model except that this assumes a circular stimulus space (e.g. orientation and direction) and therefore fixes the position of the Gaussian pRF to be at a constant distance from the origin.
6. *Reverse Correlation*: This procedure uses a model-free approach to estimate the pRF profile. It convolves the stimulus apertures with the HRF, then regresses the activity measured in each fMRI volume against these apertures, so that for each pixel in the aperture you see how much the stimulus predicts the measured response. This gives an image that reflects at which location the voxel showed the strongest responses. Because our stimuli are usually temporally ordered and predictable designs there will be a lot of correlation throughout the image (e.g. with a bar stimulus you will see a set of bars intersecting at the pRF location). Some filtering may therefore be required to minimise these artifacts whilst retaining the (presumably) true pRF. These pRF profiles are then saved and the x and y coordinates of the pRF peak and its spread are estimated from that. Up to this point the procedure is comparably fast. However, you can also fit 2D pRF models to those reverse correlation profiles. This requires you to define some additional parameters in the Model. You can also transform pRF profiles (rotating them into a common axis and/or centring them on their peak) and average them for further analysis. This procedure has many advantages over model-based pRF fits but also substantial disadvantages. So the choice between forward-model and reverse-correlation methods really depends on the purpose of your research.

7. *Connective Field Reverse Correlation*: This method is vaguely based on the ideas first presented by Haak et al. 2013 in *NeuroImage*. It estimates the vertex-wise functional connectivity (correlation between time courses) between all vertices in your region of interest, and each vertex in a pre-defined *seed region*. In a retinotopic paradigm, this seed region could for example be V1. This produces a correlation profile across the whole seed region. The connective field is the region where you observe the strongest correlation. We then use a template map of the seed region to assign a position to each connective field. This template could be the probabilistic retinotopic map predictions described by Benson et al. 2012 in *Current Biology* or some anatomical definition. Further, we can then fit a 2D pRF model to the correlation profile projected back into stimulus space to estimate pRFs.
8. *Connective Field Forward-Model*: Like for pRF models, you can fit a forward model for connective fields. This is similar to the method first described by Haak et al. 2013. This has some advantages over using reverse correlations, but it also comes with other caveats and depends on several assumptions. It also contains a few options to simplify the analysis, such as analysing only polar angle wedges or eccentricity bands.

## Surface projection

One last thing you should know about how the analysis in SamSrf works is surface projection. FreeSurfer's automatic cortical reconstruction creates a 3D model of the boundary between grey and white matter. This is a mesh consisting of over 100,000 little points (vertices) that are connected to one another in an even greater number of tiny little triangles (faces). FreeSurfer also expands this surface until it hits the pial boundary between grey matter and cerebrospinal fluid. Thus, for each vertex on the grey-white mesh there is a corresponding vertex in the pial mesh. We can use this to find points that lie within the grey matter and locate which voxels those correspond to in our functional MRI scans. This is what we do when we project the fMRI data onto the cortical surface. All analyses in SamSrf take place on the surface (but functions for running volumetric analyses are available if needed. These have never been tested by us however, and you will lose a lot of functionality).

There are different options for conducting the surface projection. You can either do that in SamSrf using the included functions, or you can do it directly in FreeSurfer. The latter approach is faster and you can include volumetric or surface-based smoothing if desired and use a number of other features. We recommend using FreeSurfer because of its speed – however, we probably cannot answer any questions about this approach, since we didn't write that code.

FreeSurfer further inflates the cortical meshes to smooth out the cortical folding structure. This is helpful for displaying the data on the cortical sheet and delineating maps. Finally, it also expands these inflated meshes to a spherical mesh. This can be used to combine data from different individuals. It is also extremely useful for surface-based (geodesic) operations that we might want to do later (e.g. smoothing, field sign calculation, cortical magnification, etc.).

## SamSrf for pRF models

This cookbook is a generic manual indicating how to conduct a simple population receptive field (pRF) experiment using SPM, Freesurfer and SamSrf for data analysis. Currently, SamSrf stands for ***Seriously Annoying Matlab Surfer***. If you come up with a better acronym that meets Sam's approval, you will win a prize and the acronym will be changed in your honour.

Note that this guide has been written in part from the perspective of a UNIX-based system user (e.g. Linux, MacOS). If you are working in Windows, working through a virtual machine may be necessary although as of Windows X it is possible to run the UNIX terminal scripts directly in Windows. This guide assumes basic knowledge of MATLAB, SPM, and of the UNIX command line.

**Important note to users of SamSrf versions prior to 6: A lot of things have changed in SamSrf 6. In particular, the GUI has been removed.** Most people will use scripting sooner or later for most of their analysis pipeline and the GUI started to become increasingly complex and it was not worth the effort supporting it any longer. There are now a few GUI tools available for displaying maps and delineating ROIs, and a documentation GUI for the parameters you specify in the Model structure.

Required dependencies:

- **MATLAB 64-bit** (*version R2020a or later  
earlier versions & 32-bit versions have not been tested*)
- **Statistics** and **Optimization** toolboxes from MATLAB (*only for some functions*)

Not required but we use:

- **SPM8** or **SPM12** *for pre-processing (but you could use any other package)  
& if using SamSrf for surface projection  
There is also a patch for loading NII files  
using MATLAB's native functions instead.*
- **FreeSurfer 6.0.0** or later – *not strictly required but it's what we use for surface models*

**Note that SamSrf saves data MAT files in format 7.3.** This should not be a problem for most people but **if your Matlab cannot read the MAT files** generated by SamSrf you need to change the default setting in **File->Preferences->General->MAT-Files**

### Overview of pRF pipeline

Before delving into SamSrf analysis proper, let's look at the overall plan for data analysis in pRF modelling. The analysis takes place in three component streams: pre-processing in an fMRI package of your choice (e.g. SPM, FSL), surface reconstruction in FreeSurfer, and pRF modelling in SamSrf. Between these streams, the whole analysis can be broken down into 5 main steps:

1. fMRI pre-processing (in any package)
2. Anatomical reconstruction (FreeSurfer)
3. Surface-projection of functional data (either SamSrf or FreeSurfer)
4. pRF modelling with inputs from all previous steps (SamSrf)
5. Displaying results, area delineation, and quantitative analysis (FreeSurfer or SamSrf)



## CHAPTER 1: ACQUIRING DATA

In order to perform a pRF experiment you will need to acquire the following MRI scans:

### Echo-planar imaging (EPI)

- One or more runs of a visual mapping stimulus
  - In this example, we use 2 runs of geometric bars in the following configuration:
    1. Bars traversing the visual field in four different directions
    2. The same bars but sweeping in a different order
- (Optional) You may collect one run of sparse visual stimulation for estimating hemodynamic response functions (HRF) but we will normally use the canonical HRF
- (Optional) You may want to collect B0 field maps to correct distortion in your EPI scans (probably more important for 3T and 7T data than 1.5T)

### Structural images

- A high-resolution T1-weighted anatomical image with the full 32-channel head-coil on.
- (Optional) If acquiring EPI scans with the top part of the head-coil off, a short T1-weighted anatomical image also with the top coil section off. Generally, when reusing an older anatomical scan for your participant, it is advisable to acquire a fast structural to aid coregistration with the old scan on which your FreeSurfer reconstruction is based.

## CHAPTER 2: DIRECTORY STRUCTURE

For each of your subjects, you will need to create an organised directory structure to perform the analysis. See Figure 1:

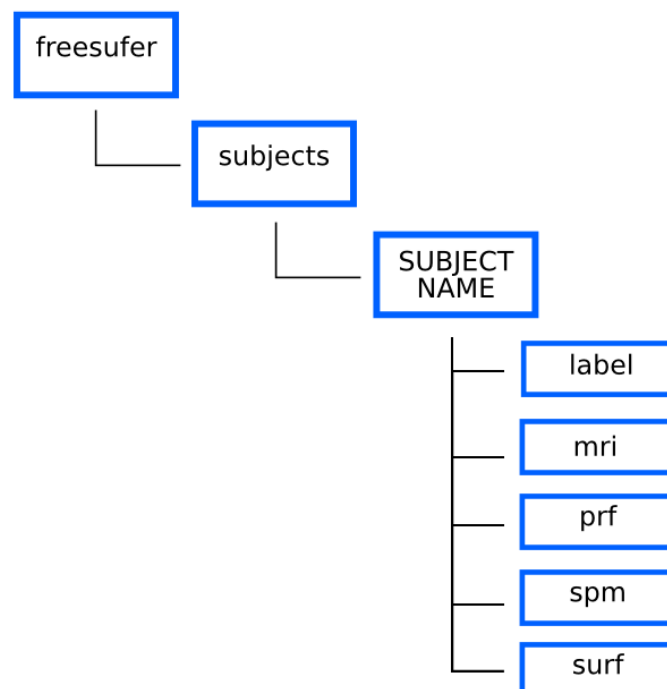


Figure 1

The analysis for each subject is contained within a directory specific to that subject. Within the subject directory you will find:

*label* contains the surface labels that FreeSurfer creates during reconstruction/parcellation  
*mri* contains the volume information and the original T1-weighted image reconstructed  
*surf* contains the surface reconstruction information from FreeSurfer  
*spm* contains the files for the pre-processing steps will reside (you can choose this name)  
*prf* is the population receptive field analysis directory (you can choose this name)  
*anatomy* is the folder where the anatomical meshes for the subject will be stored

**Note that most of your directory structure will be created when you run recon-all in FreeSurfer (step 3), except the folders *spm*, *anatomy*, and *prf*. The folders *spm* and *prf* will need to be created manually and you can choose the names for these.** (In theory, you can also keep those folders somewhere else if you so desire, but this cookbook assumes you leave them with the FreeSurfer data). The folder *anatomy* is created when you surface-project data.

If you wish to carry out more than one pRF experiment in a single subject (e.g. two conditions), then a sensible step might be to create two folders (*prf\_1* and *prf\_2*) instead of one, and treat them independently from chapter 4 onwards. Alternatively, you can also keep data from different conditions in the same folder but name the analyses with different, descriptive names.

**Note that the subjects folder shown here is inside the FreeSurfer installation. However, you can use any arbitrary folder as your *subjects* folder**, e.g. when you have different projects and you don't want the data to be mixed. To change the *subjects* folder you need edit *SetUpFreeSurfer.csh*, which is located inside your main FreeSurfer folder. You will find a line that says something like:

```
setenv SUBJECTS_DIR /usr/local/freesurfer/subjects
```

You can change the path to wherever you want to keep your data. However, you have to make sure that various folders FreeSurfer wants (from the original *freesurfer/subjects* folder) are present.

## CHAPTER 3: PREPROCESSING

### 1. fMRI pre-processing

Copy your data to SUBJECT/spm and perform your standard pre-processing for fMRI images in your software of choice.

In the case of data acquired at the FIL, we recommend the following steps in SPM8:

- DICOM import\*
- Bias correction (only when using 32-channel head coil)
- Remove dummies
- Optional: Create B0 field maps
- Realign & unwarp (with or without B0 field maps)
- Optional: Slice-timing correction (Karl Friston will tell you it's useless)
- Co-registration to the structural
- Optional: Smoothing (but this is best done on surface at cortical projection stage)
- 4D merge (in some centres this is the normal output from the scanner)

An expanded explanation of these steps is provided in Appendix A.

From your pre-processing, you should end up with a single 4D file for each stimulation run in NiFTI format. In this example, we have 3 files: 1<sup>st</sup> run of fMRI (bars\_1), 2<sup>nd</sup> run of fMRI (bars\_2) and the HRF scan.

You also need to convert the T1-weighted structural from .hdr/.img to .nii.

## 2. Run recon-all

Locate the T1-weighted file and run the automatic reconstruction:

```
recon-all -i /path/to/T1.nii -s SUBJECT -all
```

Note SUBJECT is going to be the subject name in the *freesurfer/subjects* directory.

Additionally you can use the `-label_v1` option in this command to derive an anatomical prediction of V1 (based on Hinds et al., 2008, NeuroImage). This does not work on all versions of FreeSurfer. We tested SamSrf extensively on 6.0.0 but we always recommend using the latest version.

## 3. Inspect surface

Check that the reconstruction worked in each hemisphere with *tksurfer*. Example:

```
tksurfer SUBJECT lh inflated
```

Instead you can also do this in SamSrf, using the *SurfaceProjection* tool and by projecting data to the surface and opening that in *DisplayMaps*.

## 4. Create occipital labels

This is optional but advised for saving time. There are two ways to create the occipital ROI.

1. The default in SamSrf is now to do this automatically in MATLAB using the *MakeOccRoi* command (in *SamSrf/Utils*). This reads in the surface data created by *recon-all* and restricts the occipital ROI label to the posterior vertices.
2. Originally, we created these labels manually in *tksurfer*. For the sake of completeness we leave these instructions in: Open the surface file as in step 3 and draw a wide region encompassing the occipital lobe.

Save label with meaningful names, e.g. 'lh\_occ' and 'rh\_occ' in the *prf* directory

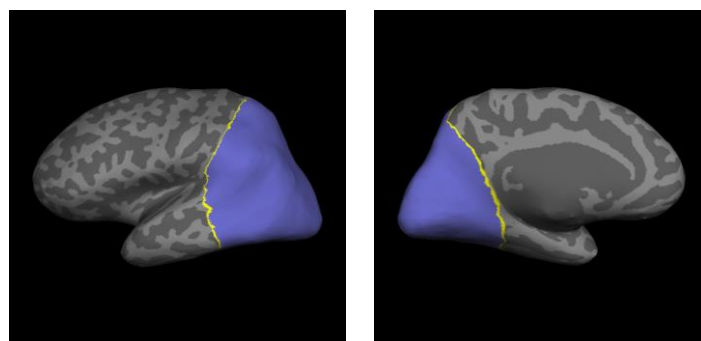


Figure 2

## CHAPTER 4: SAMSRF PROCESSING

As of SamSrf 8.2, you can also read in GII (GIfTI) files directly. You can also read in surface data projected using FreeSurfer with `mri_vol2surf`. If you are using this functionality you can select the MGH/GII files you created for each functional run using `mri_vol2surf` in `SurfaceProjection`. If you do this you probably already know how to use FreeSurfer functions. In fact, we generally use this procedure ourselves because it is much faster than our functions. Nevertheless, for completeness we continue below with the instructions using SamSrf's own surface projection functionality **but we recommend using other approaches**.

**Important:** for structural data with resolution better than 1mm (and some 7T data in general) the SamSrf functions currently do not work – so you must use FreeSurfer's functions or another way of surface projection instead!

### 5. Check structural registration

In MATLAB, navigate to where you want the pRF data to live (usually in `SUBJECT_NAME/prf`)

Now launch `SurfaceProjection`


Select the main folder for your FreeSurfer subject (labelled SUBJECT NAME in Figure 1)

After a moment of loading the files a SPM check registration window will open

A red line indicating the grey-white matter boundary will appear in both hemispheres

*Note: since SamSrf 6.30 you can also use Matlab's native NII reading functionality by using the `NatMatlabNii4SamSrf` patch. This cannot display the surfaces interactively like SPM so you will have to make do with what it can do. Only use this if you know what you're doing.*

Look at the match-up between the anatomy and the red overlay as shown in Figure 3. It should define a boundary between the white and grey matter that follows the sulcal pattern of the cortex:

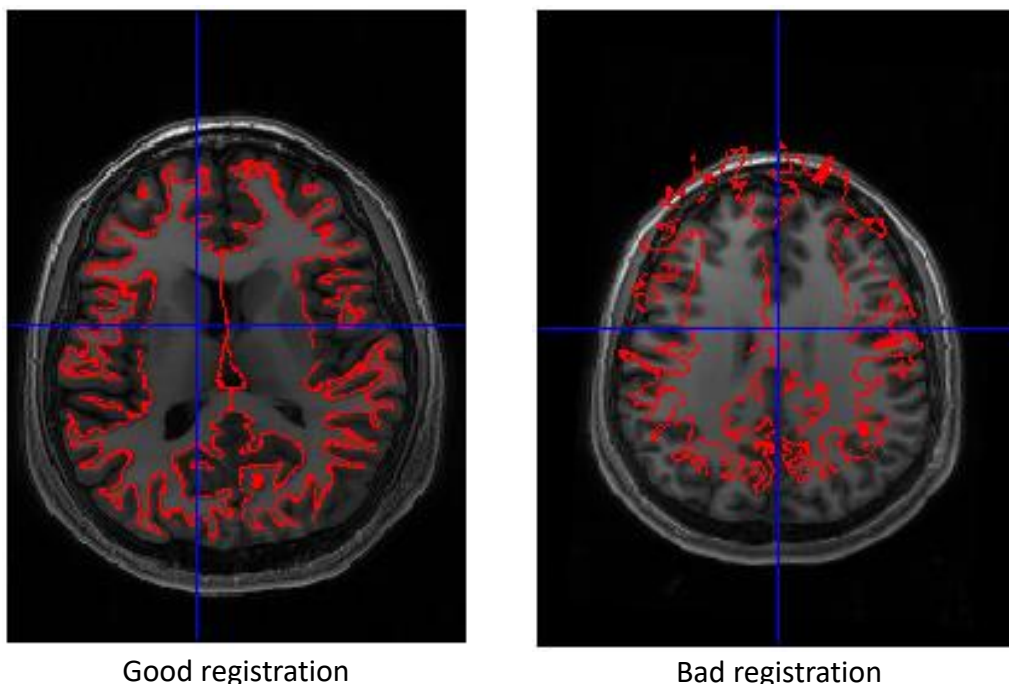

Good registration

Bad registration

Figure 3A

There are various reasons why the link between FreeSurfer's coordinate system and SPM may fail. For example, if you fiddled with the affine transformation of the structural scan you used for FreeSurfer reconstruction (or if you resliced it or used it as source image in coregistration), then the alignment will be broken.

The example in Figure 3A should be pretty obvious but note that things can be a lot more subtle, as in Figure 3B. Note that the surface is only off by a few millimetres:

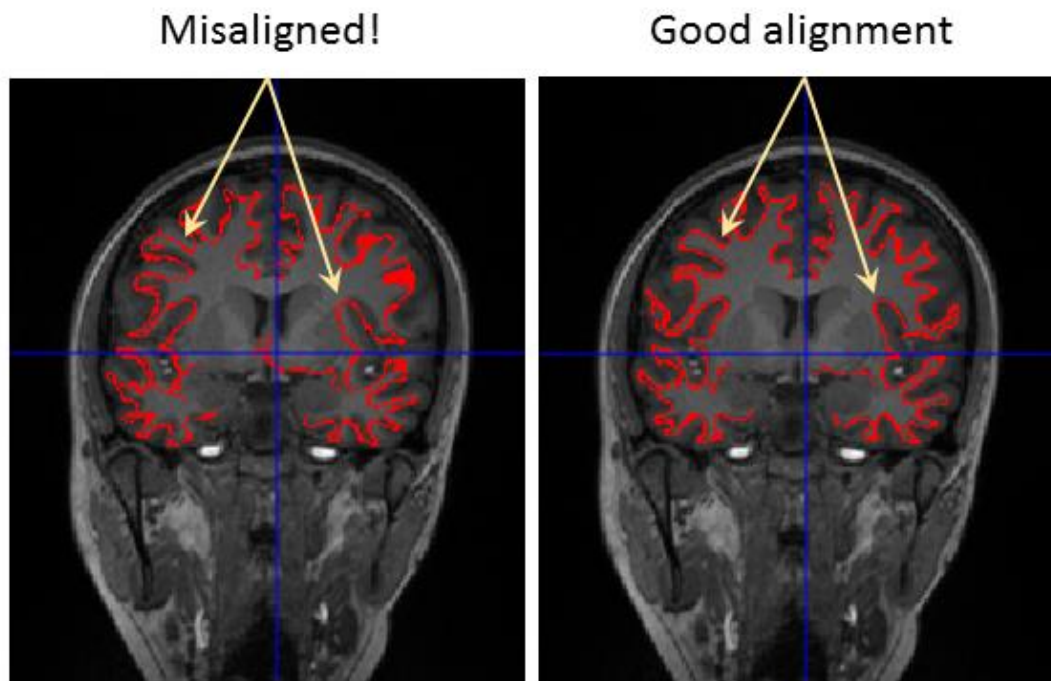


Figure 3B

This tool is also helpful for inspecting the surface reconstruction has worked properly (instead of step 3 above). If the red surface doesn't cover the whole cortex but parts of the cortex are outside of it, then you will have to fix the reconstruction. This may be due to too sharply defined structural images. Subtly smoothing (~1mm kernel) the structural in SPM can help in this situation.

Note that if you use `mri_vol2surf` then step 5 is only to ensure the surfaces are correct.

.....

### What to do if registration is bad

In previous versions of SamSrf, there was a hack where you could provide a coregistration file that ensured the alignment between FreeSurfer's and SPM's space coordinates works. As of SamSrf version 6.3 this has been removed because we believe to have fixed it. **Generally, for T1 data with better than 1mm resolution and/or ultra-high field data you should really use FreeSurfer's functions – and we recommend using those in any case!**

.....

## 6. Create surface data files

SurfaceProjection will then automatically ask you for the other required files:

Select FreeSurfer subject folder (e.g. *freesurfer/subjects/SUBJECT/*)

This will look for the .nii file from freesurfer that was used for reconstruction in *freesurfer/subjects/SUBJECT/mri/NAME\_OF\_T1.nii* or its orig subfolder.

Select functionals: choose all the pre-processed 4D NiFTI images, e.g. bars\_1.nii, bars\_2.nii  
Note that if you use mri\_vol2surf then select the MGH files for your functional runs instead.

You can also provide the *Cortical sampling steps* as a vector. The number is in proportions of cortical thickness. If you have more than one value here, then the activity of different steps will be averaged. If you don't enter anything, a default of 0.5 is used. Note that if you use mri\_vol2surf then this step is not available as you have already done the cortical sampling.

Next there is a dialog asking you whether you want to *average*, *concatenate* or *separate* runs. Obviously, for averaging the stimulus protocol must then be identical in all those runs. For concatenation you may run into problems if there are too many runs to put into one matrix which will depend on the amount of RAM in your machine. Finally, you can keep individual runs separate (the projection will be done separately for each then and thus take much longer). For most purposes, averaging or concatenating are more useful but we left this function in for completeness' sake.

Next there is a dialog asking whether you want to use temporal normalisation. For functional data analysis, normalisation is necessary but you may also want to project all sorts of other data (e.g. structural data, GLM betas, t-statistics, etc) so then you will not want to use it.

Finally, there is also a menu allowing you to specify the projection method. This menu only appears if you defined more than one cortical sampling step. By default, this is the mean across cortical sampling steps but it can also be the median, minimum, maximum, or you can turn it off so all sampling steps are saved (this is useful for plotting the data along the vertex normal – **currently the main analysis functions in SamSrf however do not support such data**).

This will take a few minutes (or longer if you have many cortical sampling steps). Once it's done, SamSrf will output all the files you gave it in .mat format with the prefix 'lh\_' or 'rh\_' appended to indicate which hemisphere it belongs to. The SurfaceProjection tool does this automatically for both hemispheres.

**Remember to do this for both the functional runs from pRF mapping as well as the HRF measurement run (if you have one).** The functions used to extract the HRF also work in surface space so they need this data. **You can also combine surface projected data of both hemispheres into a bilateral (bi\_\*.mat) file using the samsrf\_bilat\_srf.m function.**

*Note:* as of SamSrf 6.1 the surface projection stage also calculates the noise ceiling of your data – but only when using averaged runs. This is based on splitting the data into odd and even runs and so works best when having an even number of runs. The noise ceiling gives

you an upper bound of what goodness-of-fit your pRF model can possibly achieve with these data. It is stored in `Srf.Noise_Ceiling` in raw data files but after model fitting it is moved into `Srf.Data`.

## 7. Fit HRF

**This is optional: you can simply use a canonical HRF or fit the HRF using your pRF data**

In MATLAB, run this function to fit an HRF to the data from your HRF run. Assuming your HRF run is called `lh_uabfHrf.mat` and the TR of that run was 2.55 seconds you run:

```
samsrf_extract_hrf('lh_uabfHrf', 2.55)
```

Check the curve-fit is sensible, for example (Figure 4):

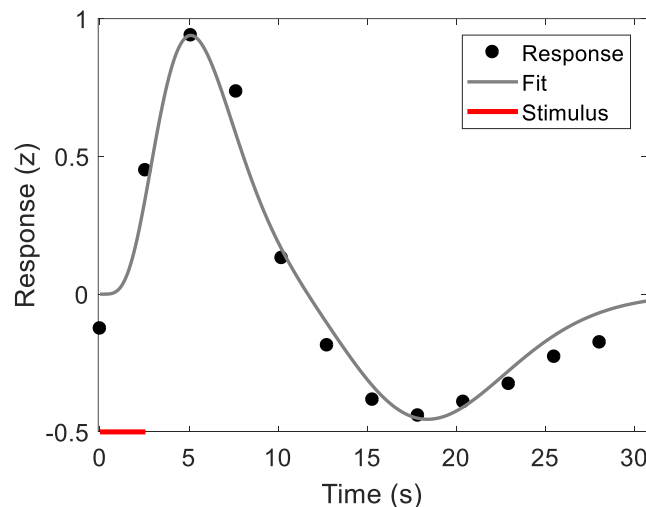


Figure 4

*Fix bad HRF:*

This step is rarely ever necessary. If your HRF data contain outliers that skew the HRF fit, you can fix this using the `samsrf_fix_broken_hrf` function.

## 8. Create apertures

Next you need the apertures to tell the model where the stimulus was in each volume.

If you used a standard retinotopy design using some of the SamPenDu lab's scripts, you should be able to do this using the `CreateApertures` function and load the results files containing the apertures. If you used your own stimulus protocol, you will have to generate them in your own way.

Critically, the apertures should be stored in a .mat file in a variable called `ApFrm`.

This is a 3D matrix where columns and rows correspond to the x and y coordinates of pixels, and the third dimension is time. There is a frame for every volume. Usually these are binary images (1 for stimulus, 0 for background) but that is not essential.

**Important: For forward-modelling, you must vectorise these apertures using `VectoriseApertures`, which will save them with the `_vec` suffix!**

## 9. Inspect apertures

Once generated/saved, you may want to check that the apertures look as they should.

Run the `ViewApertures` tool.

Select the aperture file (should be named something like *aps\_\*.mat*).

Click through the display to make sure everything is correct.

You can also click *Play Movie* to play it automatically.

## CHAPTER 5: pRF MODEL ANALYSIS

Unlike in previous versions, since SamSrf 6 there is no longer a GUI for running the analysis. More people scripted their analysis pipeline anyway and the GUI wasn't really that convenient. We now have wrapper scripts that define the parameters of the analysis. You can adapt them for your batch scripts to completely automate the process.

Example scripts are kept in the *SamSrf/Models* folder. There are scripts for the basic pRF models you might want to use, such as the standard 2D Gaussian, the Difference-of-Gaussians, various asymmetric pRFs, as well as 1D tuning curve models. You can also create your own – but this is way outside the scope of this cookbook. Here we will describe what you need to do to run a standard 2D Gaussian model defined in the script *Standard\_2d\_Gaussian\_pRF.m*

The idea is that you take a copy of these example scripts, put that in a place related to your analysis, such as the parent folder where you keep all your subject data or a folder where you keep your other analysis scripts. You then modify the script to define the parameters of your analysis. It might also be advisable to give your script a unique name that reflects the project. For instance, in a study on pRFs in aging we might want to call this function *Aging\_Study\_2d\_pRF.m*. In the script, you define a structure called *Model* which contains various parameters. If that seems scary, don't fret because most of these you won't need to worry about unless you define your own pRF model. You can simply skip forward to step **15-18** and then on to **20**.

**Importantly, you can use the `ModelHelp` tool to learn about all the possible model parameters in the various analyses and how to use them.**

## 10. Define the pRF model

*Model.Prf\_Function* is a MATLAB functional handle to the pRF model function. For the standard 2D Gaussian pRF model this is defined like this:

```
@(P,ApWidth) prf_gaussian_rf(P(1), P(2), P(3), ApWidth)
```

The handle has two inputs: *P*, which contains the actual model parameters to be fitted, and *ApWidth*, which is an internally required variable defining the stimulus coordinates.

The maximum number of parameters is ten. You can see how other models are defined in the other example scripts in *SamSrf/Models*. **You don't need to worry about this part unless you want to define your own pRF models.**

## 11. Model name

*Model.Name* defines the name of the pRF model, which is what the data file will be called. In our example script, this is 'pRF\_Gaussian' so the final map file for the left hemisphere



will be called *lh\_pRF\_Gaussian.mat*. You can change this if you want.

## 12. Parameter names

*Model.Param\_Names* is a cell array defining the names of the parameters. In a retinotopy model, the first two will be the horizontal and vertical coordinate of the pRF centre, so we call them *x0* and *y0*. For the pRF size (standard deviation of the Gaussian), we define the name as *Sigma*. The standard 2D Gaussian model only has three parameters so we can stop there. If we fit more parameters, they need to be named here. **The beta (amplitude) and intercept (baseline) are fit separately, after the pRF parameters have been estimated.**

## 13. Determine scaled parameters

Depending on the receptive field model you use, some or all the parameters are relative to the eccentricity of your mapping stimulus. Some parameters (like *x* and *y* coordinates or pRF size) are in this space and must therefore be scaled. The model fitting will reject parameter estimates twice the *Model.Scaling\_Factor* (see below).

*Model.Scaled\_Param* is a Boolean vector and it toggles whether a parameter is scaled or not. In the case of the standard 2D Gaussian model, all three parameters are scaled so it is [1 1 1]. But in the case of some other models, such as the oriented asymmetric 2D pRF, one of the parameters is the orientation of the pRF and this is obviously *not* scaled. **Critically, this vector must be of equal length as the cell array in step 12.**

## 14. Determine positive-only parameters

The model fitting function will look at the fitted (optimised) parameters in turn and reject those with absurd estimates. Any parameters that are scaled (see 13. *Model.Scaled\_Param*) will automatically be rejected if the estimate is greater than twice the scaling factor. However, some parameters should also only be positive, e.g. the pRF size or tuning width cannot be negative.

*Model.Only\_Positive* is a Boolean vector that toggles whether a parameter must be positive. In the case of the standard 2D Gaussian model, this is only the pRF size (*Sigma*).

## 15. Define scaling factor

Once we know which parameters will be scaled, we also need to define by how much it will be scaled. In two-dimensional retinotopic models, the scaling factor is the maximal eccentricity of the mapping stimulus (usually the radius of the mapped area). In one-dimensional tuning curve models (for example, tonotopy) this would be the half-width of the stimulus space because the middle is defined as 0.

**You must define this critical parameter because this will inevitably differ for your setup!**

## 16. Define the TR

You will also need to define the TR of your pulse sequence in seconds. This is critical because the convolution of the predicted time series with the HRF depends on the TR.

**This is also a critical parameter that depends on your experimental setup!**

## 17. Specify the HRF

In *Model.Hrf* you will also need to specify which HRF to use. We typically use the canonical HRF in which case we leave this empty []. If you estimated the HRF as described earlier in step 7, you can provide here the file name of the estimated HRF. Or you can provide the

HRF directly as a vector, where each component corresponds to one TR. Finally, if you don't want any HRF to be used (as you might in some situations) you must set this to 1. **Whatever you do, you need to specify this, even if it just [] to use the canonical.**

## 18. Specify aperture file

Finally, in *Model.Aperture\_File* you will need to specify the file containing the apertures. Since this will typically depend on your experimental setup, **you must define this**. If you used the exact same stimulus design for each participant, you could just keep the aperture file in a common folder (with your model script) and provide the full path name here. **Apertures must be vectorised and contain both the ApFrm and ApXY variables!** You can use *VectoriseApertures* to convert movie apertures into vectorised ones.

## 19. Define search space

You can also define the parameters of the search space used for the coarse fitting (extensive grid search) stage. For each of the five parameters you define a vector of the points of the search grid for that parameter, so this vector determines both the range and the granularity of the search space. **Unlike in previous versions, since SamSrf 9 the search grid must be defined in stimulus space (e.g. degrees of visual angle)!** You could define it in unit space and multiple the scaled parameters by the *Model.Scaling\_Factor* (see example scripts in SamSrf/Models)

For most people, there probably isn't much of a reason to change these parameters. If you make the granularity too fine the search space will become unwieldy and you may even run out of memory. The following parameters can be defined. These are just names for the model engine. The actual names in the final map you already defined above in step 12. The parameters can mean several different things, depending on the pRF model. In models with fewer than ten parameters you can just omit unneeded ones (they are internally set to 0).

<i>Model.Param1:</i>	The horizontal coordinates of 2D pRF centres or the polar angle location of 2D pRF centres or the peak location of a 1D tuning curve mode
<i>Model.Param2:</i>	The vertical coordinates of 2D pRF centres or the eccentricity of 2D pRF centres or the tuning width of a 1D tuning curve model
<i>Model.Param3:</i>	pRF size, usually the standard deviation of a Gaussian pRF
<i>Model.Param4:</i>	Second pRF size in the DoG pRF or elliptical pRFs
<i>Model.Param5:</i>	The rare fifth parameter. In the DoG model the amplitude ratio, in oriented model this is the orientation, etc...
<i>Model.Param6-10:</i>	More complex models will have even more parameters although it becomes questionable how accurate these model fits will be. Ten parameters is probably more than the algorithm can realistically cope with...

## 20. Run model fit

Now we are ready to run the model fitting. This is done using the generic *samsrf\_fit\_prf* function which takes the following inputs: the Model structure you specified above, a cell array with the SamSrf data files we created in Chapter 4, and (optionally) the region of interest label we created in step 4 (see the examples in SamSrf/Models). More complex pRF models or 1D tuning curves are run in a similar way. You can hopefully deduce a lot of that

from the example scripts and by referring to the `ModelHelp` tool.

## OPTIONAL MODEL PARAMETERS

In addition to these standard parameters there are quite a few optional parameters, and analyses like reverse correlation of connective field modelling require different parameters. Please refer to the `ModelHelp` tool for more information about these.

### 21. Post-processing pRF map

The pRF fitting function does not use any smoothed data unless you tell it to smooth the coarse fit or you presmooth your data before running the model fit.

You may also wish to consider denoising the map. Beta parameters can be very diagnostic about rubbish data. If the Beta is 0 or negative, or if it is massively large, this is a sign of useless data. Using the `samsrf_denoisemap` function you can set the  $R^2$  of such vertices to 0 and thus filter them out. Similar filtering may also be necessary based on implausible Sigma estimates (for example really tiny pRFs below the resolution of your stimulus or which are not biologically plausible).

Further, for delineation or even just making figures you may wish to have smoother maps. You may also want to calculate field sign maps or the cortical magnification factor. To do all of that we can load the pRF map file and run the function `samsrf_surfcals`. In addition, for more complex models you may need to do other post-processing, for example calculating the FWHM of the DoG model or ordering orientation parameters. This can also be done at this point and some of the example scripts show how you might do such post-processing.

## RECAP: WHAT FILES DO I NEED TO RUN MODEL FITTING?

Here is you need to have in the *prf* folder (change *lh* to *rh* for right hemisphere):

- Apertures file** with same length as the whole time series: (e.g. *aps\_Bars.mat*)
- Surface data files for **pRF mapping runs** (e.g. *lh\_Bars1.mat*, *lh\_Bars2.mat*, etc)
- (Optional) HRF parameter file from an **HRF run** created in step 7 (e.g. *hrf\_lh\_Hrf.mat*)
- (Optional) **ROI label** file to restrict analysis created in step 4 (e.h. *lh\_occ.label*)

## INTERMISSION: WHAT HAPPENS DURING MODEL FITTING?

First the surface data files are loaded into memory. If a region-of-interest was chosen the mask is loaded. Next, the search space is generated (see the *Introduction to pRF models* section). This used to take a long time on earlier machines but nowadays this is probably fine for most people. For more complex models (e.g. the Difference-of-Gaussians pRF) this takes longer because the search space is more complex. If you have the exact same order of conditions for each subject you can reuse the same search space (*src\_\*.mat* files) for each though and thus save time.

Next, the search space predictions are convolved with the haemodynamic response function. Then the coarse fitting starts. On modern computers and MATLAB versions this should take no more than a few minutes. Then the fine fitting commences and this can take from a half hour to a day or even a week, depending on what model you're fitting, which procedure you chose, and which

MATLAB version you are running. Finally, the maps are saved in a .mat file. **Note that the Model variable in which you defined the model parameters (such as the search space dimensions and the HRF) is also saved in that file! Be aware of this if you are reloading the .mat file!**

## CHAPTER 6: RENDERING MAPS

### 22. Displaying maps

To inspect maps and save figures you can use the DisplayMaps tool. This will ask you for a pRF map file in SamSrf .mat format and then load this up. By default this limits the view to an occipital ROI (you can change this option in *SamSrf\_defaults.mat*).

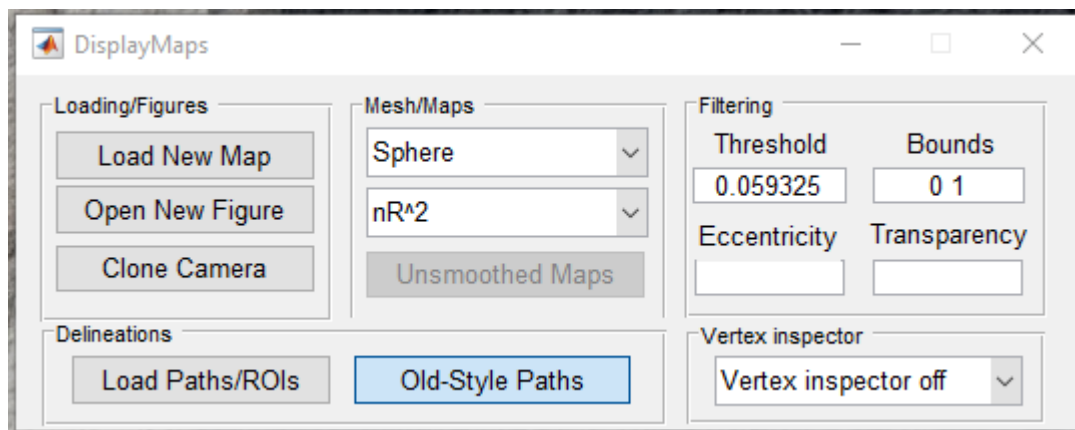


Figure 5

You can load any of the maps stored in your SamSrf file. The Polar Angle and Eccentricity maps need to be calculated as the file only stores Cartesian X and Y variables.

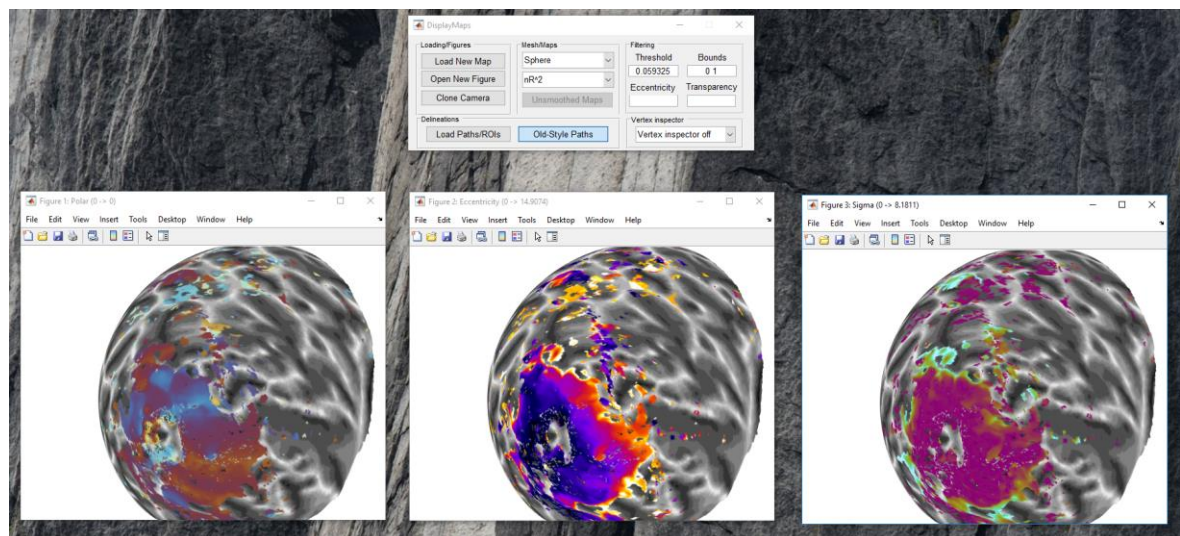


Figure 6

You can also load labels or paths previously saved. These can be in FreeSurfer .label format (ASCII files with a list of vertices) or paths in .mat files. You can also load paths stored in a delineation file saved by the DelineationTool

You can select any of the following cortical meshes to display: *Inflated*, *Pial*, *White*, and

*Sphere* (plus some other options derived from these). The spherical mesh is the default because it allows you to see a lot more cortex in the same view, but it may not be ideal for knowing what you're looking at (especially outside of occipital cortex).

The *Threshold* is automatically adjusted to give you a pre-defined p-value for your  $R^2$ . You can also manually adjust this to your liking. Similarly, you can adjust the *Eccentricity* range as a two-element vector (e.g. 1 8 to only show eccentricities between 1° and 8°). The *Bounds* determines how the maps are clipped. For  $R^2$  maps this value determines at which level the hottest/colds colours are displayed. For eccentricity maps, this value determines beyond which eccentricity the maximum/minimum colour is used etc.

Finally, you can also set the *Transparency*. This determines the transparency of the map. By default this is 0.1, which means that that all model fits above the 10<sup>th</sup> percentile are fully opaque but below that the transparency is ramped down to the threshold. You can also set this to a negative value. In that case, the whole map has the same transparency, so for instance -0.5 means that the map is 50% transparent.

You can use the *Open New Figure* button to create a new main figure. This is the only figure that is updated by anything you do in the tool. So, you should first load the maps you want and then open a new figure for any further maps you want to look at simultaneously. If you have several figures open (e.g. the polar, eccentricity, sigma and  $R^2$  maps) and want them all to show the same view, you click the *Clone Camera* button. This means all the other plots will adopt the same view as your main figure.

You can load a completely new map using the *Load New Map* button.

## 23. Renderer lighting

The pretty lighting can be nice for figures but does not seem to work on all computers and may result in a display bug. On some computers, the display is also very slow when the lighting is turned on. You can turn the lighting in the figures on and off with this command:

```
samsrf_lighting('off')  
samsrf_lighting('on')
```

Lighting is turned off by default. For pretty figures, you may want to **first** adjust the camera angle and then turn on lighting.


## 24. How to delineate regions

A detailed description of how to delineate visual regions (including the higher extra-striate cortex) is the subject of a separate cookbook. In short, visual regions are bounded by the reversals in the polar and eccentricity maps. For the early regions (V1-V3) the polar map is fully sufficient. For higher regions (e.g. V4 and in particular the V3A/B complex and VO-1/2) it helps to also look at the eccentricity map because most of these regions have foveal representations separate from the early visual areas. A field sign map (showing whether or not a map is a mirror image of visual space) can also be useful, although the clarity of this map can vary a lot. Finally, the beta (amplitude) map can be very helpful in identifying bad data, in particular off-the-edge artifacts in parts of the visual field that were not stimulated. The *AutoDelineation* tool allows you to create some automatically drawn borders for regions

based on an atlas. These are far from perfect and must be corrected/completed, but they improve the efficiency and reproducibility of delineation.

To help you with the delineation Figure 7 shows the colour wheels. They can also be found in the *SamSrf/Colours* folder in your SamSrf installation or you can display them in MATLAB using the function

```
samsrf_colourcode
```

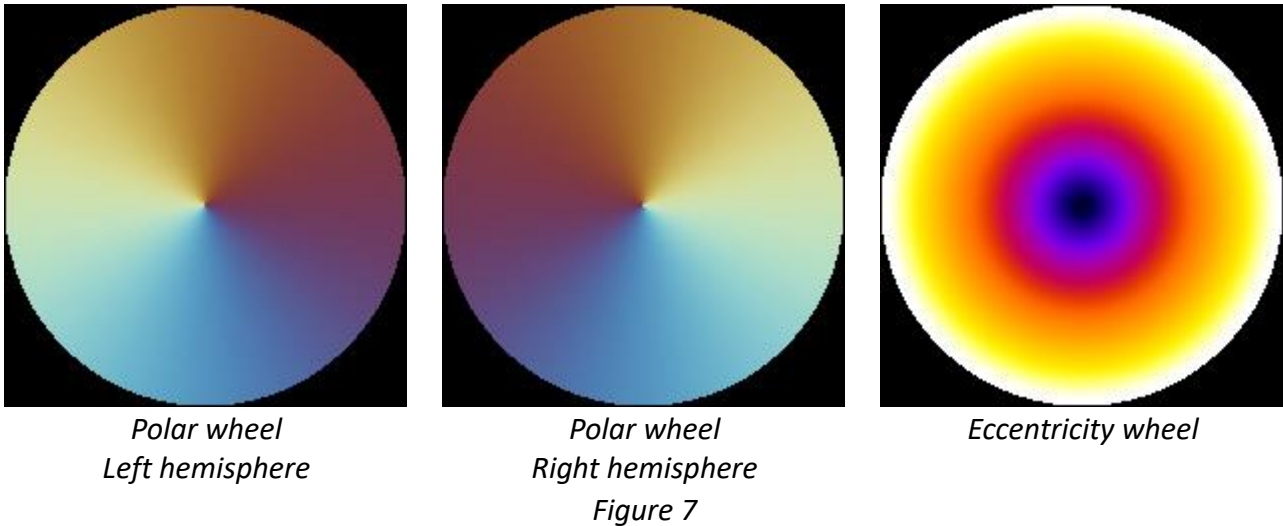

## 25. Exploring the results

To save the variables that generate the above graphs, you need to run these command line functions and repeat for each ROI:

```
load lh_pRF_Gaussian  
Srf = samsrf_expand_srf(Srf)
```

This loads the pRF data into the workspace as a variable Srf. **The second line is critical to expand the data set and add in the anatomical meshes! Without that, none of the other functions may run properly**, although all native functions in theory should do that automatically (see also Appendix B).

To then plot data for a particular region of interest you can run:

```
samsrf_plot(Srf, 'Sigma', Srf, 'Eccentricity', [0 9], 'lh_V3', 0.5)
```

This produces a scatter plot for sigma in Srf against the eccentricity in Srf (Figure 8). The vector [0 9] determines the eccentricity range to be plotted. You can provide two different Srf variables here, e.g. when you want to compare different conditions based on a common eccentricity estimate. Here 'lh\_V3' is the name of a ROI label for V3 in the left hemisphere. The 0.5 refers to the  $r^2$  threshold to be used ( $r^2$  is  $R^2$  normalised by noise ceiling).

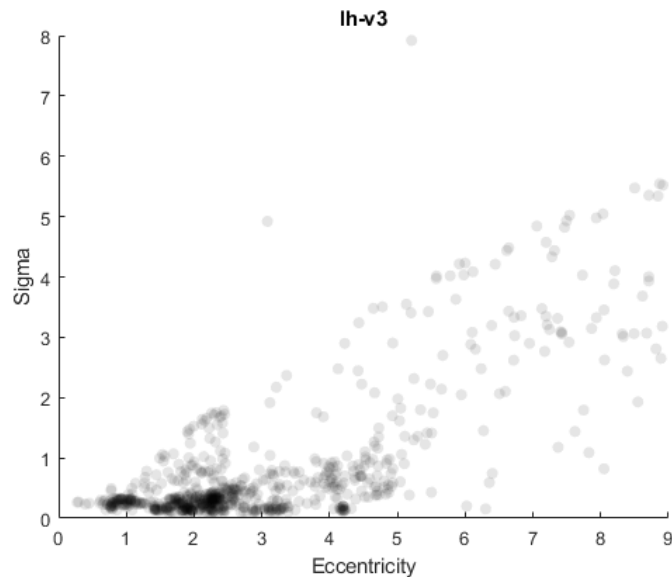


Figure 8

You can also use binned plots including polar plots when using polar angles for binning. **But please be aware of the dangers of binning analyses!**  
 \_\_\_\_\_(see <https://www.biorxiv.org/content/10.1101/2020.12.15.422942v2>)

There are many other functions available for analysing and quantifying the pRF results. For instance, you may want to plot the visual field coverage of V1 for pRFs with the  $r^2$  threshold  $> 0.5$ , assuming a mapping stimulus with eccentricity of 9, and clipped to  $>10\%$  of pRF coverage (Figure 9).

```
samsrf_vfcoverage(Srf, 9, 'lh_v1', 0.5, 0.1)
```

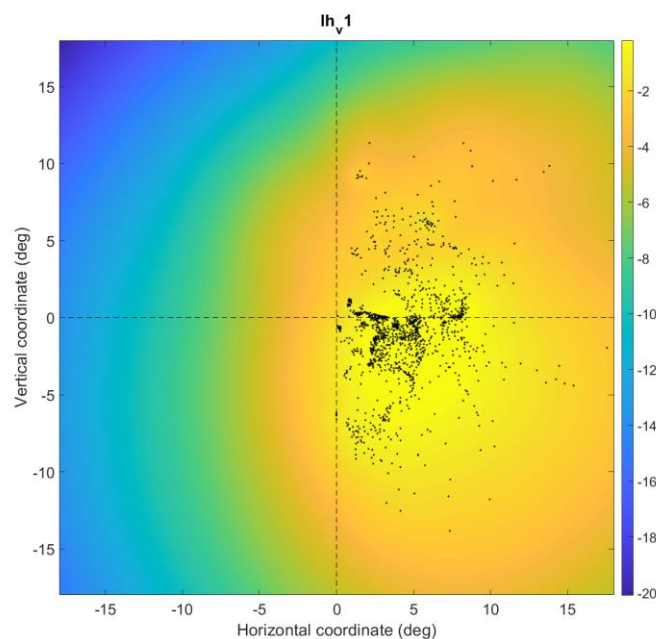


Figure 9

All functions contain help information that should explain how they work.

## APPENDIX A: FMRI PRE-PROCESSING

Sample pre-processing steps for fMRI data carried out in SPM8. **Note that these instructions are specific to data acquired in the FIL on a 3T TIM Trio scanner with a 32-channel coil.** If you didn't acquire your data there under these conditions, the procedure for you will most likely differ. A lot has changed since these notes were written, but they are still correct in principle.

### DICOM import

SPM > DICOM import

Convert all your raw files to 3D NiFTI format, conserving the original filenames

*prefix added: [f]*

### Bias correction

Perform bias correction on all the EPI images using spm\_biascorrect

*prefix added: [b]*

### Remove dummies

Place the first four volumes in a separate folder

### Create fieldmap vdm\_\* files

After DICOM import, fielfmaps are divided into two folders; the first has two files (magnitude) and the second has a single file (phase).

Batch script

SPM > Tools > Fieldmap > Presubstracted phase and magnitude data

Phase image: second field map folder (single file)

Magnitude image: first file in first field map folder

Blip direction: -1

EPI-based field map?: no

Jacobian modulation?: no

Brain mask = no

Match to the first volume in each EPI run

### Realign and unwarp

Batch script

SPM > Spatial > Realign & Unwarp

1 subject = 1 module, each run in a separate session

for each run, use the single session-matching vdm\_ file

Masking = yes

*prefix added: [u]*

### Slice-timing correction

Batch script

SPM > temporal > Slice timing

On all bias-corrected, realigned and unwarped functionals

Insert scan details (TR, TA, number of slices, slice order)

Reference slice is middle one (15 of 30)

*prefix added: [a]*

If you acquired a single structural, do a one-step coregistration. If you acquired a short and a long



structural, do a two-step coregistration. Note that it is important to select only 'estimate' and not 'estimate and reslice' in order to preserve the spatial origin from the fMRI images throughout.

### One-step coregistration

Batch script – Estimate

All functionals to high-resolution structural

Reference (target): high-resolution structural from Freesurfer folder

Source: 1<sup>st</sup> aubf\_\* image of 1<sup>st</sup> run

Other images: all other aubf\_\* volumes for all runs

One module per participant

### Two-step coregistration

Batch script

SPM>Spatial>Coreg>Estimate

All functionals to short structural with coil off

Reference (target): short structural

Source: 1<sup>st</sup> aubf\_\* image of 1<sup>st</sup> run

Other images: all other aubf\_\* volumes for all runs

SPM>Spatial>Coreg>Estimate

All functionals to high-resolution structural with coil on

Reference (target): high-resolution structural from Freesurfer folder

Source: short structural

Other images: dependency on results from step above

One module per participant

### Check coregistration

SPM>CheckReg

### Optional: Smoothing (use only with the pre-smoothed receptive field model!)

All the coregistered, slice-timing corrected, unwarped, and bias-corrected volumes

Kernel size: 4mm

*prefix added: [s]*

### 4D merge

Combine smoothed images for each run using spm\_file\_merge

Good idea to rename them here into the correct run, for example:

- bars\_1
- bars\_2
- HRF

## APPENDIX B: STRUCTURE OF SAMSRF DATA FILES

### Introduction

You carefully designed your study, pre-processed your data and ran a model fit with SamSrf. Now, you can use the buttons in `DisplayMaps` to plot maps or used some functions to plot data. But what if you want to get your hands on the data themselves?

This guide covers the contents of the SamSrf results files and their potential uses. It assumes basic knowledge of the MATLAB environment and programming language. If you wish to learn MATLAB, you can head to MathWorks<sup>1</sup> or the excellent guide of MATLAB for Psychologists<sup>2</sup>.

### Navigation

To access the SamSrf results file, open MATLAB and go into the subject folder where you ran the pRF fit. Here you will find a file called *lh\_pRF\_Gaussian.mat* or similar.

- 'lh' or 'rh' for left and right hemisphere, respectively, or 'bi' for bilateral files, or 'vol' for a volumetric data file (feature largely untested!)
- Or any arbitrary name for the actual data prior to any analysis (e.g. 'lh\_Bars1')

Open the file by double clicking or typing:

```
load('lh_pRF_Gaussian.mat')
```

Once loaded, you will see 2 variables in workspace (some formats may contain more):

Model:	The pRF model specified as described in Chapter 5
Srf:	The main structure with all the data

By default, Srf variables are compressed and have their anatomical meshes removed. So before you can use them you must expand them (but all of the native SamSrf functions should in theory do that automatically). To do so run:

```
[Srf, roi] = samsrf_expand_srf(Srf)
```

Finally, before saving .mat files with a Srf you may wish to compress them again:

```
Srf = samsrf_compress_srf(Srf, roi)
```

This will remove the anatomical meshes from the Srf. Further, it also restricts the Srf to a region of interest defined by the vector *roi* (e.g. the occipital region). See *samsrf\_compress\_srf* for details.

---

<sup>1</sup> [http://www.mathworks.co.uk/academia/student\\_center/tutorials/launchpad.html](http://www.mathworks.co.uk/academia/student_center/tutorials/launchpad.html)

<sup>2</sup> <http://antoniahamilton.com/matlab.html>

## **Contents of 'Srf'**

The variable Srf is a structure containing several sub-variables. Note that this list only describes the most basic fields you might see when running a forward-model pRF analysis. Other analyses will contain additional fields:

**Srf.Version:** [scalar] Contains the version number used by the analysis. As of SamSrf 6, this is first added during surface projection but then updated again during pRF model fitting. Also, this has only been introduced in version 4.0 so if an analysed map file doesn't have this field it was analysed with an older version. This field is important to ensure data are comparable (e.g. Sigmas before version 4 were measured incorrectly). But ideally you shouldn't have to use any files created by versions prior to 7...

**Srf.Structural:** [string] Full path location of the structural image used as reference in alignment and for cortical reconstruction in FreeSurfer.

**Srf.Functional:** [string or cell array] If this is a pRF map, this contains a description of analysis carried out (e.g. which pRF model) and details of smoothing (if any) and name of ROI restricting the analysis (if any). If this is a data file, this contains the name of the functional NII file the data came from. If this is an average, it is a cell array of all the file names that were averaged.

**Srf.Hemisphere:** [string] Hemisphere contained in this file. This should match the filename. Please note that you can also create bilateral fields containing both hemispheres, in which case this will be 'bi'. You can also run volumetric analyses in which case this will be 'vol'.

**Srf.Cortex\_Steps:** [scalar] Cortical sampling steps, as defined in the SamSrf GUI. This defines the voxel positions and number that each vertex is sampled from in volumetric space and is in units of cortical thickness. Therefore, the default 0.5 would mean that each vertex contains the data from that voxel that is located 50% along the normal vector from the vertex on the grey-white matter surface and the corresponding vertex on the pial surface. You could also choose more voxels though by setting this to a vector, e.g. [0:.2:1] would sample six voxels along this vector between grey-white and pial surface in steps of 20% of the cortical thickness. Taking more than one voxel probably only makes sense for relatively high-resolution data. Note that when you use FreeSurfer's surface projection functions this field is uninformative.

**Srf.Meshes:** [char] Contains a link to refer to the folder and file name where the anatomical mesh data are stored. This is not mandatory but is used by default when analysing data since SamSrf 6. It means that when you have more than one file with functional data projected to the surface, you only need to have one file for the anatomical meshes rather than duplicating this across files. Obviously, if you only have one file (e.g. for a public data repository) you might just share that with all the necessary data in the same file.

**Srf.Vertices:** [ $m \times 3$  matrix] An index list of every vertex. Each row corresponds to the X, Y and Z coordinates of a single point in native space using FreeSurfer's internal coordinate system. The vertex coordinates correspond to the grey-white matter surface, so gyri and sulci are folded. May be stored in a separate file for anatomical meshes (see *Srf.Meshes*).

**Srf.Faces:** [ $m \times 3$  matrix] Face values for every vertex. The surface meshes are constructed from thousands of tiny triangles, called *faces*. Each face is defined by three vertices and they are defined by the vertex indices in this list, one face per row. May be stored in a separate file for anatomical meshes (see *Srf.Meshes*).

**Srf.Normals:** [ $m \times 3$  matrix] Each rows contains the vector pointing from a vertex on the grey-white matter surface to the corresponding vertex on the pial surface. This is used by surface projection to determine the voxels to sample from. May be stored in a separate file for anatomical meshes (see *Srf.Meshes*).

**Srf.Inflated/Sphere/Pial:** [ $m \times 3$  matrix] Like *Srf.Vertices* but contains the coordinates model for the inflated, sphere, or pial surface mesh instead of the white-grey matter boundary. May be stored in a separate file for anatomical meshes (see *Srf.Meshes*).

**Srf.Curvature/Area/Thickness:** [ $1 \times m$  vector] Contains the curvature, surface area, or thickness of the cortical surface model for each vertex. May be stored in a separate file for anatomical meshes (see *Srf.Meshes*).

**Srf.Data:** [ $k \times m$  matrix] This contains all the resulting values of the analysis. If this field *Srf.Values* (see below) is present, the rows correspond to the  $k$  value descriptions in *Srf.Values* (e.g.  $R^2$ ,  $X_0$ ,  $Y_0$ , Sigma, etc..) and each column corresponds to a single vertex. So, if you are interested in, say, the  $R^2$  values of every single vertex, you can index the first row and every column of *Srf.Data*. If *Srf.Values* is not present this is probably a file containing the data prior to analysis and each row of *Srf.Data* contains one volume. If you applied surface smoothing post-analysis, this variable contains the smoothed data. The un-smoothed data is stored in *Srf.Raw\_Data*.

**Srf.Rule:** [char] Defines how the voxels in *Srf.Voxels* relate to the data for each vertex in *Srf.Data*. This only makes a difference if you have more than one cortical sampling step. The default is 'Mean' which means that each vertex is the arithmetic mean of the sampled voxels. Other options are 'Sum' for the sum, or 'Median' for the median (see *samsrf\_vol2srf* for details). This is obviously only relevant for the actual surface projection. The analysed data files inherit this field as this tells you how surface projection was done.

**Srf.Noise\_Ceiling:** [ $1 \times m$  vector] Contains the noise ceiling if it has been calculated. This is only the case for raw data files. If the noise ceiling exists, the model fitting procedure will add this to *Srf.Data* and so there won't be a separate field in *Srf* for that in those files. At the moment, surface projection functions (*samsrf\_vol2srf* and *samsrf\_mgh2srf*) will automatically calculate the noise ceiling if there is an even number of runs that are being averaged. This is because the noise ceiling is based on split-half reliability data. When you concatenate runs or have more complicated designs, you may need to derive the noise ceiling yourself.

**Srf.Y:** [ $n \times m$  matrix] The preprocessed time series for each vertex after z-normalizing each run and concatenating them (or in the case of CF analysis, after global mean correction). Each column is for a vertex, each row is for a volume (time point). This only exists for analysed data files.

**Srf.X:** [ $n \times m$  matrix] The *neural* (rather than metabolic, so it is not convolved with the HRF!) time series for each vertex predicted by the pRF model. Therefore,  $R^2$  can be calculated by comparing *Srf.X* to *Srf.Y*, after convolving *Srf.X* with the HRF. **Important: Just to complicate things, when running only a coarse fit, this field contains convolved time series!**

**Srf.Values:** [ $k \times 1$  cell array] Description of the values stored in *Srf.Data*. Each row in this variable corresponds to the matching row in *Srf.Data*, so if the 4th row is 'Sigma', all the values along the 4th row in *Srf.Data* are Sigma values. Used for indexing. As of SamSrf 6, this is mandatory (and automatically added by *samsrf\_expand\_srf*, if it doesn't exist).

**Srf.Raw\_Data:** [ $p \times m$  matrix] If you applied surface smoothing, this variable will contain the non-smoothed data. If no smoothing is applied, this variable will not exist. Because values like the CMF are calculated after smoothing, they are not contained in this field.

**Srf.Roi:** [vector] A list of vertex indices that indicate whether any given vertex is inside the ROI used to limit the analysis. If used, the number of columns in *Srf.Data* should match the length of *Srf.Roi*. If this is the case, use *samsrf\_expand\_srf* to expand the data. This is essential before running any further analysis scripts on the data.